


AirLife™ Closed Suction System

Critical Care Ventilation

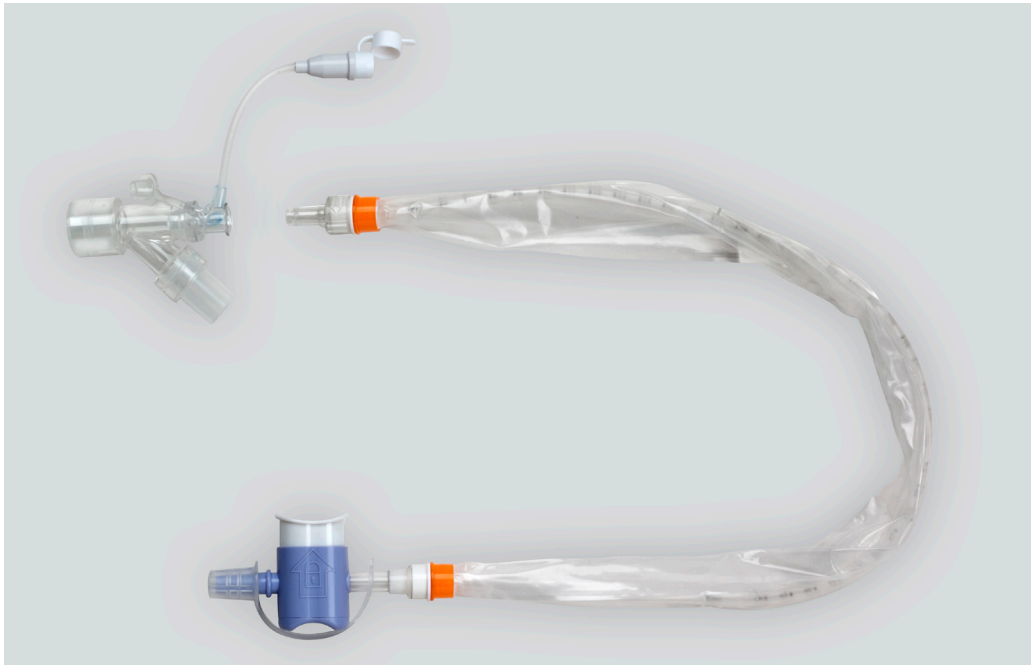

Promote *safety*
and demand
versatility

Introducing the Verso™ Airway Access Adapter

The most versatile adapter
for safe and simple
access to your ventilated
patient's airway



Raise the bar, reduce the risk

Clinicians understand that disconnecting or “breaking” of the ventilator circuit can lead to lung de-recruitment and hypoxemia, and it creates possible opportunities for cross-contamination. The AirLife™ Closed Suction System is a multi-purpose device made up of modular components, helping to keep your patients safer than traditional closed-suction catheters. The system features the Verso™ Airway Access Adapter and its revolutionary Pucker™ Valve technology that seals the patient's airway from the access port. The Verso™ adapter is placed in-line when setting up the ventilator circuit, so you can perform closed suction catheter change-outs and various airway access procedures through a single port **without ever disconnecting the circuit.**

Safety: Pucker™ Valve helps maintain safe, continuous ventilation while you change closed suction catheters and perform various other airway access procedures

Versatility: The Verso™ adapter gives you multi-purpose airway access to perform:

- Open suction
- Closed suction
- Bronchoscopic procedures
- Mini-BAL
- Instilled drug delivery

Protect all of your ventilated patients with one closed system solution

Adult/Pediatric (Verso™ Adapter) and Neonatal (Neo-Verso™ and Neo-Verso™ ‘Y’ Adapters) options provide a full suite of products with closed suction catheter sizes ranging from 5 french to 16 french and specialty catheters including tracheostomy-length, directional-tip and fluid/access catheters.

AirLife™ Closed Suction System

Living, breathing innovation to improve ventilated patient outcomes

Adult/Pediatric products*

Description	Order Number	Quantity
Adapters		
Verso™ Adult/Pediatric Airway Access Adapter	CSC100	20/box
Closed Suction Catheters for use with Verso™ Airway Access Adapters		
10 French Catheter	CSC110	50/case
12 French Catheter	CSC112	50/case
14 French Catheter	CSC114	50/case
16 French Catheter	CSC116	50/case
Specialty Catheters for use with Verso™ Airway Access Adapters		
14 French Directional-Tip Catheter	CSC114D	50/case
10 French Tracheostomy-Length Catheter	CSC110T	50/case
12 French Tracheostomy-Length Catheter	CSC112T	50/case
14 French Tracheostomy-Length Catheter	CSC114T	50/case
16 French Tracheostomy-Length Catheter	CSC116T	50/case
Adapter Seals for Bronchoscopes for use with Verso™ Airway Access Adapters		
Adapter Seal for 3.8mm-4.4mm sized Bronchoscopes	CSC100B2	25/box
Adapter Seal for 4.8mm-5.5mm sized Bronchoscopes	CSC100B3	25/box
Adapter Seal for 5.9mm-6.3mm sized Bronchoscopes	CSC100B4	25/box

Neonatal/Infant products*

Description	Order Number	Quantity
Adapters		
Neo-Verso™ Neonatal/Infant Airway Access Adapter	CSC200	20/box
Neo-Verso™ 'Y' Neonatal/Infant Airway Access Adapter	CSC300	20/box
Closed Suction Catheters for use with Neo-Verso™ Airway Access Adapters		
5 French Catheter	CSC205	50/case
6 French Catheter	CSC206	50/case
8 French Catheter	CSC208	50/case
Specialty Catheters for use with Neo-Verso™ Airway Access Adapters		
5 French Fluid/Access Catheter	CSC205S	50/case
Adapter Seals for Bronchoscopes for use with Neo-Verso™ Airway Access Adapters		
Adapter Seal for 2.8mm-3.4mm sized Bronchoscopes	CSC100B1	25/box


*All system components of the AirLife™ closed suction system may not be commercially available as of the publication of this brochure.



Adult AirLife™ Closed Suction System



Closed suction control valve

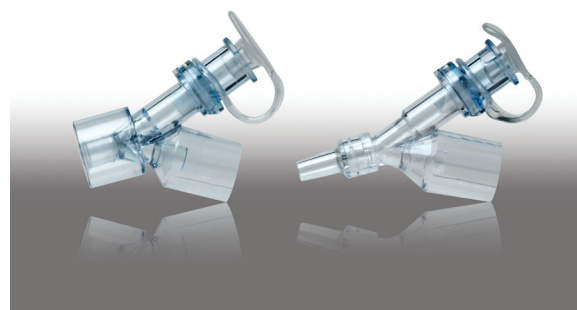

Neo-Verso™ and Neo-Verso™ 'Y' adapters

⚠ WARNING

- Single patient use only
- Do not sterilize or re-sterilize
- Federal (U.S.) law restricts this device to sale by or on the order of a physician or licensed practitioner



Manufactured for:
CareFusion
Yorba Linda, CA



Authorized Representative:
CareFusion Germany 234 GmbH
Leibnizstrasse 7
Hoechberg 97204, Germany

CareFusion
Yorba Linda, CA

© 2009 CareFusion Corporation or one of its affiliates. All rights reserved. AirLife, Neo-Verso, Pucker and Verso are trademarks of CareFusion Corporation or one of its affiliates. RC0246

